

Gel Etch

Conveniently packaged in two 5 cc syringes, Gel Etch provides an easily applied alternative to traditional etching procedures. Disposable syringe tips are included to provide a contamination-controlled environment between patients.

Gel Etch Kit 740-0150


Contains:

Syringes (two 5 cc)

Disposable Syringe Tips (25)

Syringe (one 5 cc) 740-0133

Disposable Syringe Tips (25) 740-0134

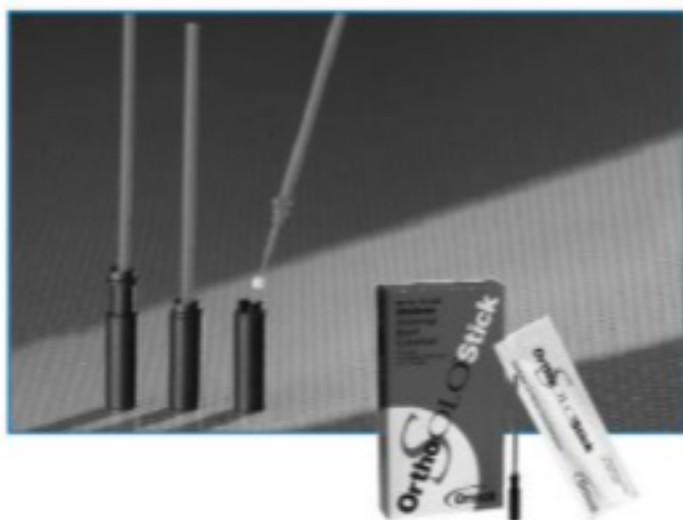

Ortho Solo™ Stick

Your Favorite Bond Enhancer in a Convenient Unidose™ Delivery System

Each Ortho Solo Stick contains enough sealant/bond enhancer for one arch, making everyday bondings faster and easier than ever before. Simply push the Ortho Solo Stick into the sealant well to access the bond-enhancing solution, pull the stick out and immediately apply the liquid to the tooth. No air-dry or light-cure steps needed.

Consistent and accurate dosage minimizes product waste. Unidose delivery reduces the risk of cross contamination and offers a cleaner presentation to your patients. A tightly sealed container ensures that there is no harmful exposure to light and no risk of product evaporation.

Box of 10 Sticks (1.2 ml total)

740-0269

Ortho Solo

Sealant and Bond Enhancer

Ortho Solo™ is a fluoride-releasing universal sealant and bond enhancer. With light-cure and paste-paste systems, Ortho Solo replaces the sealant. With no-mix chemical-cure systems, Ortho Solo should be used in addition to the activator for best results.

Ortho Solo incorporates a bond-enhancing property that improves adhesion to the tooth at the adhesive interface, hence reducing bond failures. The glass filler, unique to Ortho Solo, acts as a stress and shock absorber, preventing cracks that can lead to bond failures. This formulation also releases fluoride where you need it the most at the enamel-resin interface.

With Ortho Solo, no air-dry or light-cure step is needed, thereby saving valuable chair time. It also includes a compound that displaces moisture, making the technique more forgiving when a dry field is inadvertently compromised.

Ortho Solo Kit 740-0270

Contains:

Ortho Solo (four 5 ml bottles)

Applicator Brushes (75) 740-0041

Disposable Wells (1 pkg)

Ortho Solo Technique Card 070-0440

Ortho Solo Refill 740-0271

Contains:

Ortho Solo (one 5 ml bottle)

Liquid Etch

Etching Solution

Contains 37% phosphoric acid.

Bottle of 10 ml

740-0038