

## Ordering Information

ANTERIORS	Torq.	Ang.	D/O	Part Number
◆ Upper 3-3 Universal	+55°	0°	0°	369-2102
◆ Upper 3-3 Universal, w/ hook	+55°	0°	0°	369-3102
◆ Lower 3-3 Universal	+40°	0°	0°	369-2110
◆ Lower 3-3 Universal, w/ hook	+40°	0°	0°	369-3110

BICUSPIDS	Torq.	Ang.	D/O	Part Number
Upper/Lower 4 & 5 Universal	+11°	0°	0°	369-2400
Upper/Lower 4 & 5 Universal, w/ hook	+11°	0°	0°	369-3400
Upper/Lower 4 & 5 Universal	0°	0°	0°	369-2401
Upper/Lower 4 & 5 Universal, w/ hook	0°	0°	0°	369-3401

MOLARS	Torq.	Ang.	D/O	Part Number
Upper 6R, with hook	+10°	0°	+5°	369-3600
Upper 6L, with hook	+10°	0°	+5°	369-3601
Upper 7R, with hook	+10°	0°	+5°	369-3700
Upper 7L, with hook	+10°	0°	+5°	369-3701
Lower 6R, with hook	0°	0°	0°	369-3610
Lower 6L, with hook	0°	0°	0°	369-3611
Lower 7R, with hook	0°	0°	0°	369-3710
Lower 7L, with hook	0°	0°	0°	369-3711

PATIENT EDUCATION/MARKETING	Part Number
STb Lingual Typodont	717-1040
STb Social 6 Patient Brochures, English (40 per pack)	070-5493
STb Social 6 Patient Brochures, Spanish (40 per pack)	070-5494
STb Social 6 POP Display, English (brochure holder)	070-5495
STb Social 6 POP Display, Spanish (brochure holder)	070-5496
STb Social 6 Poster, English	070-5497
STb Social 6 Poster, Spanish	070-5498
STb Demonstration DVD	070-5555

◆ Parts for STb Social 6 cases.

SINGLE PATIENT KITS	Part Number
Extraction Upper & Lower 7-7	740-0411
Upper & Lower 7-7	740-0412
Upper Only 7-7	740-0414
Upper Only 7-7	740-0415
◆ Upper 3-3 Only	740-0416
◆ Upper & Lower 3-3	740-0417

ARCHWIRES – PACKS OF 10	Part Number
<b>Leveling 1</b>	
◆ .013 Copper Ni-Ti (8.5" Straight Length)	266-1360
◆ .013 Preformed Copper Ni-Ti (STb Social 6)	205-1950
<b>Leveling 2</b>	
◆ .016 Ni-Ti-Upper Size 1	205-0025
◆ .016 Ni-Ti-Upper Size 2	205-0023
◆ .016 Ni-Ti-Upper Size 3	205-0021
◆ .016 Ni-Ti-Lower Size 1	205-0031
◆ .016 Ni-Ti-Lower Size 2	205-0029
◆ .016 Ni-Ti-Lower Size 3	205-0027
<b>Detailing</b>	
◆ .016 TMA	202-0025
.0175 x .0175 TMA	202-0018
.017 x .025 TMA	202-0020



# Ormco | AOA

# STb™

Light Lingual System



## INVISIBLE TREATMENT 3-D CONTROL OUTSTANDING RESULTS

While removable aligner treatment is attractive to many patients who demand an "invisible solution," many cases require the high level of control that only fixed appliances can provide. With the STb Light Lingual System, you no longer need to worry about patients removing and/or losing appliances, thus interrupting treatment.

The STb System combines maximum 3-D control and complete invisibility, making it the ideal aesthetic solution for you AND your patients. Exceptional quality, efficiency, simplicity and patient comfort. You can have it all with the STb System.

- Low-force, low-friction mechanics for faster tooth movement
- Maximum patient comfort with minimal impact on speech
- Simplified lab procedures
- Predictable high-quality results

Contact your Ormco sales representative to order today!

### United States and Canada

ormco.com  
Call: 800.854.1741 or 714.516.7400

### Australia and New Zealand

ormco.com  
Australia – Call: 1800 023 603  
New Zealand – Call: 0800 446 140

### Latin America

Ormco América Latina  
ormcoamericalatina.com  
Call: +52 (55) 5280 7851

Ormco de México  
ormco.com.mx

Call: +52 (55) 5280 7851

### Asia and the Middle East

ormco.com/index/ormco-company  
for your local Ormco Distributor

### Europe

ormcoeurope.com  
Call: 00800 3032 3032

## AOA

AOA Orthodontic Appliances  
P.O. Box 725  
Sturtevant, WI 53177  
aolab.com

## Ormco

Ormco Corporation  
1717 West Collins Avenue  
Orange, CA 92867  
ormco.com

© 2012 Ormco Corporation 070-5345 Rev.C

Twins Digital Auxiliaries Practice Development Education  
Self Ligation Aligners Tubes/Bands Archwires Lab Products

ormco.com | aolab.com

# A New Look at Lingual



Dr. Giuseppe Scuzzo




Dr. Kyoto Takemoto

Even for the inexperienced lingual user, the advanced technology of the STb Light Lingual System makes lingual treatment remarkably easy-to-use and more effective than ever.

Developed with Drs. Giuseppe Scuzzo and Kyoto Takemoto, the STb appliance is small, comfortable and efficient. The system is designed to utilize light forces, provide excellent 3-D control and eliminate the need for patient compliance.

The complete STb Light Lingual System is ideal for 7-7 cases, while STb Social 6™ corrects minor-to-moderate anterior crowding or spacing within a few months.

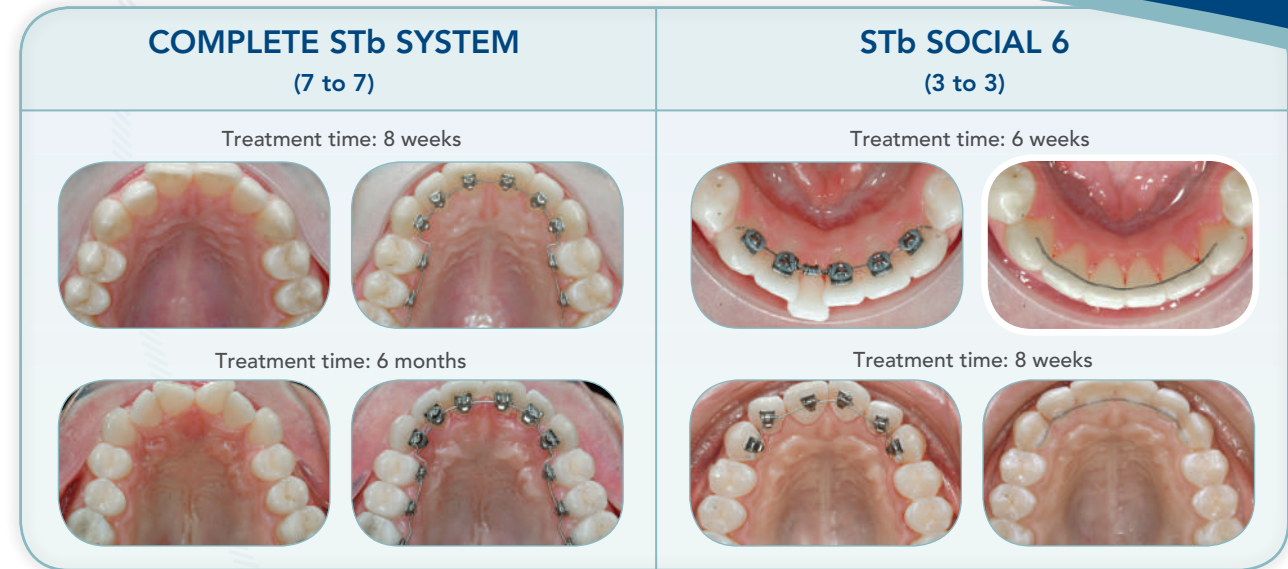

## STb Bracket Improvements

- Improved bracket ligation notch facilitates a more passive system for lower friction in larger diameter archwires
- All brackets feature a hook (option on anteriors and bicuspid)
- Anteriors – single tie-wing design at the incisal edge offers greater patient comfort; gingivally offset pad allows for more incisal bracket placement
- Bicuspid – wider tie-wing for improved rotational control and patient comfort
- First molars – now a tube for greater patient comfort



The reduced profile and rounded contours of the STb bracket minimize tongue discomfort and pronunciation challenges.

## Quality

- Treat a wide range of cases with maximum 3-D control, delivering optimal results
- Confidently predict treatment times
- Eliminate compliance issues and misplaced aligners

## Patient Comfort

- Low 1.5 mm bracket profile and rounded contours provide maximum comfort
- Patients regain proper pronunciation within just a few days
- Light-force mechanics are designed to make treatment more gentle to the tissues

## Efficiency



- Passive ligation design and optimal-force wires designed to reduce friction and accelerate treatment
- Small and narrow bracket design reduces occlusal interference and increases interbracket distance for lighter forces
- Bracket base bonds securely to tooth surfaces

## Simplicity

- Standard indirect setup in your own practice or through AOA Lab makes bracket placement fast, simple and accurate
- Most cases do not require a TARG or CLASS setup, making STb nearly half the cost of other lingual setups
- Easy for doctor and staff members to learn and use
- Available in single patient kits

Without the compliance issues and limited 3-D control of aligner treatment, the STb System delivers the aesthetics and comfort your patients will love and the clinical ease and results you demand.

It's lingual like you never seen before.

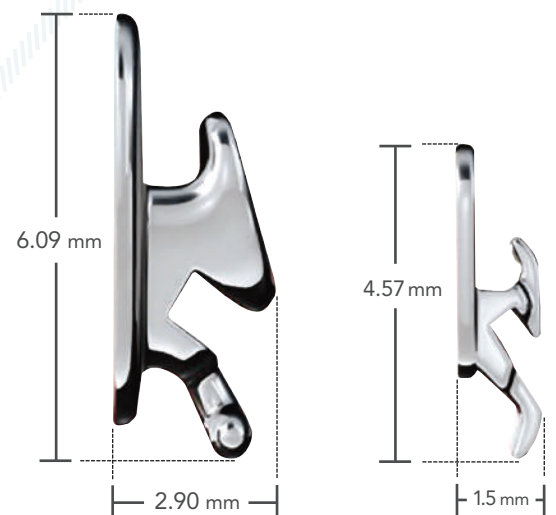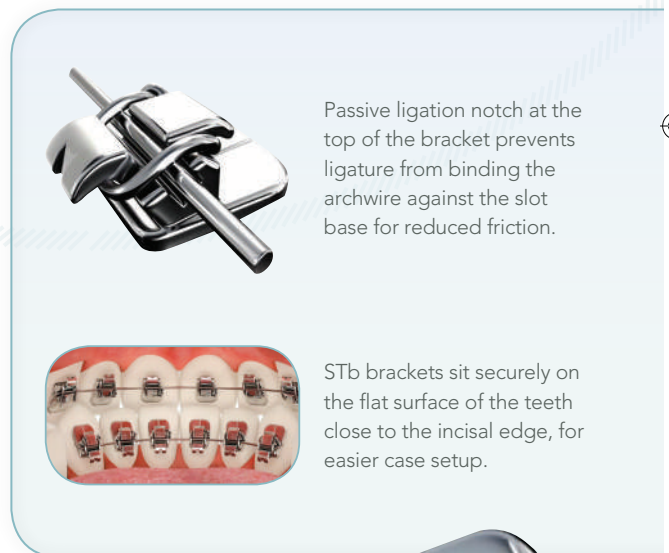

Earlier Lingual Bracket

STb Social 6 Lingual Bracket

## The Truly Invisible Aligner Alternative



- Treats minor-to-moderate anterior crowding and spacing with maximum control
- Completes most cases in 6 weeks to 4 months
- Simple two-wire sequence quickly unravels and aligns teeth; most cases finish with just one wire
- Indirect bonding through AOA lab makes bracket placement fast, simple and accurate – at a cost that's 80% less than a full case of removable aligners

